**CARRELLO CLINICO
MEDICAL TROLLEY
CHARIOT MÉDICAL
CARRO MÉDICO
MEDIZINISCHER WAGEN
CARRINHO MÉDICO
ΤΡΟΧΗΛΑΤΟ ΝΟΣΗΛΕΙΑΣ
МЕДИЦИНСКА КОЛИЧКА
ZDRAVOTNICKÝ VOZÍK**

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Gima 45718 - 45720 - 45722

Jiangsu Saikang Medical Equipment Co., Ltd
No. 35 Lehong Road, Modern Agriculture Demonstration Park,
Zhangjiagang City, Jiangsu Province, China
export@saikangmedical.com
Made in China

SUNGO Europe B.V.
Fascinatio Boulevard 522, Unit 1.7, 2909VA Capelle aan den
IJssel, The Netherlands

SK-AT75077A2, SK-ET75077A, SKR054-ET

Gima S.p.A.
Via Marconi, 1 - 20060 Gessate (MI) Italy
gima@gimaitaly.com - export@gimaitaly.com
www.gimaitaly.com

Name: Medical Trolley

Product codes:

- 1) SKR054-ET – emergency trolley
- 2) SK-AT75077A2 - anesthesia Trolley
- 3) SK-ET75077A – emergency trolley

Intend Use:

Intended to be used for help to provide immediate support in the event of health emergencies (resuscitation, defibrillation) or to contain and transport medical equipment and devices in critical contexts or contribute indirectly to the diagnosis, prevention, monitoring, treatment or mitigation of a disease or injury


CAUTIONS:


- In order to use the trolley safely, it is necessary to perform regular safety inspections on the trolley. It is recommended to check it every six months to ensure that the connections are not loose and all operations are normal.
- When the trolley frame ages and it reaches a certain servicing time, which the metal part of the trolley can be recycle.
- Avoid scratching the panel with sharp-angle appliances or knives during use, and clean it frequently to keep it clean and dry.
- If the panel is accidentally stained with stains, it is recommended to clean it. Do not use alkaline or corrosive chemicals to clean the trolley, which will cause the stainless steel surface to rust.
- Regularly check the casters to avoid hitting or overloading, which will cause damage to the casters.
- It is forbidden to force the casters to brake when braking, which will cause damage to the brake system.
- The wear of the tire surface of the wheel can be detected by visual inspection. Some thread spun yarn and other accumulations may be wrapped around the wheel, remove the bolts and nuts of the wheel, clean up the debris, and check whether the bearing of the wheel is damaged. If the parts are not damaged, can be reassembled and used.
- Judging whether the caster is replaced according to the time used, wear and tear, you can replace the wheels and make sure that the axles bolt and nuts are well put together and connected. If possible, use a lock washers or locknuts to help you with this process.
- For the casters equipped with brakes, it is necessary to regularly check whether the brakes are working properly. Check the brakes once a day or before each use. For a trolley equipped with multiple brake casters, only one brake caster can be locked at a time which will allow you to try and push the trolley also check whether one of the brake performance of each caster is good. If the function of the brake fails due to wear or damage to the wheels, immediately replace the wheels and then retest the brake system.
- If the mechanism of the brake system of the caster is damaged, and the brake needs to be repaired or replaced, please contact our after-sales service department or authorized dealer provide by our company. Every time the brakes are replaced, the brake performance of the casters must be retested.

Use method of casters

1. Lock the caster: After pushing the trolley to the desired position, press the button and lock the caster when you hear a “click”, as shown in the figure below.

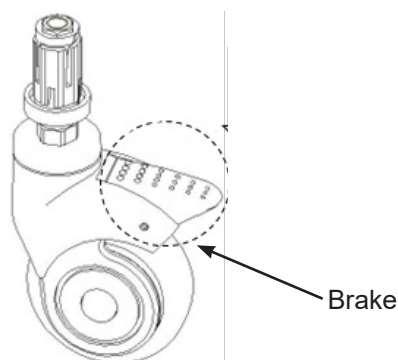2. Unlock the caster: Lift the button with your foot and unlock the castor when you hear a “click”.

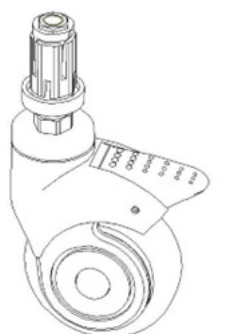Storage:

Avoid rushing out, violent vibration, and protection from sunlight and rain during the transportation of the trolley. The trolley should be stored:

- a) Ambient temperature: $-40^{\circ}\text{C} \sim 70^{\circ}\text{C}$.
- b) Relative humidity: $\leq 95\%$.
- c) Atmospheric pressure: $50 \text{ kPa} \sim 106 \text{ kPa}$.

The max load is 50Kg

After-sales service

1. Please keep the files that comes with the machine and invoices of this product properly, and you need to present these files when the company performs warranty and maintenance for the product.
2. If there is any problem in the process of using, please contact our company in time, so that our company can provide you with accurate and fast technical support and maintenance services in a timely manner.
3. From the date of sale, if the product is broken or damaged due to the correct installation and use according to the specification, the product will enjoy one-year free warranty and lifelong maintenance service with the "certificate" or invoice.
4. From the date of purchase, if it is indeed damaged or does not work normally due to quality problems within one year, the company will provide to the users a free repairing of the product.
5. Lifetime service from manufacturer: JIANGSU SAIKANG MEDICAL EQUIPMENT CO LTD

Symbol	Meaning	Symbol	Meaning
	Lot Number		Medical device compliant with Regulation (EU) 2017/745
	Warning: Read and follow the instructions (warnings) for use carefully.		Dispositivo medico
	Store away from sunlight		Manufacturing date
	Store in a cool, dry place		Unique Device Identifier
	Manufacturer		Authorized Representative in the European Union
	Read the instructions for use		Product code
	Humidity limit		Temperature limit
	Imported by		Atmospheric pressure limit

GIMA WARRANTY TERMS

The Gima 12-month standard B2B warranty applies.